


TAQ OPTIMIX CLEAR

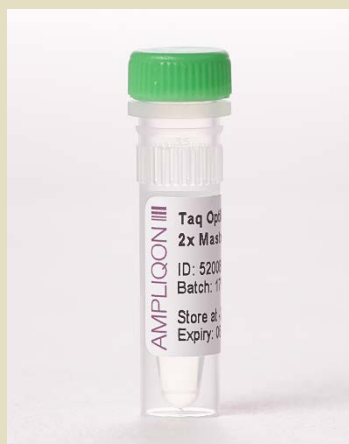


Taq OptiMix CLEAR

– An optimised Taq master mix with increased specificity

- All-in-one Taq master mix for great convenience
- Timesaving reaction setup
- Improved specificity
- Increased reproducibility

Taq OptiMix CLEAR is an optimised Taq DNA polymerase master mix, which ensures increased specificity and improved PCR performance. Taq OptiMix reduces reaction setup time and also eliminates the risk of contamination. Fewer handling steps lead to better reproducibility, making Taq OptiMix suitable for standard PCR applications, screening and high throughput testing.

Comparison with three leading Taq master mixes

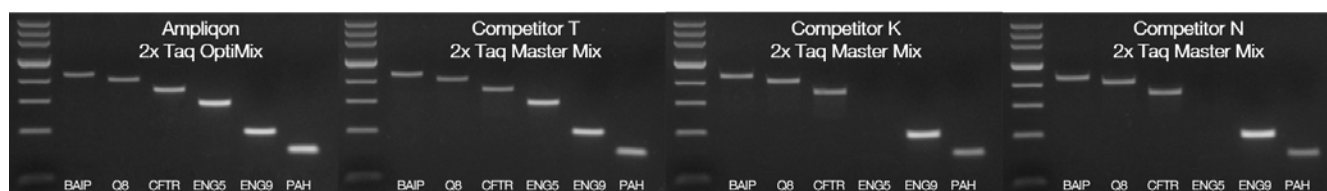

Comparison of Ampliqon Taq OptiMix CLEAR performance with equivalent Taq DNA Polymerase master mixes from three different competitors; T, K or N. Six different primer targets (BAIP, Q8, CFTR, ENG5, ENG9 and PAH) were evaluated. Taq OptiMix CLEAR performed equally well or better on all six primer targets tested in this set up.

	Size Reaction size: 50 µl	Content	Cat #
Taq OptiMix CLEAR 2x Master Mix 1.5 mM MgCl ₂ final	100 Reactions 500 Reactions 2500 Reactions 5000 Reactions	2 x 1.25 ml 10 x 1.25 ml 50 x 1.25 ml 25 x 5 ml	A370501 A370503 A370506 A370507

Taq OptiMix CLEAR 2x Master Mix

1.5 mM MgCl₂ final concentration

MADE IN DENMARK

Cat. No.: A370599 – SAMPLE

20 Reactions

	Taq OptiMix CLEAR 2x Master Mix, 1.5 mM MgCl ₂ final concentration
ID No.	5200800
Cap colour	Green
Content	1 tube x 0.5 ml

Key Features

Taq OptiMix CLEAR is a ready-to-use 2x master mix with the Ampliqon Taq DNA polymerase, the NH₄⁺ buffer system, dNTPs and magnesium chloride present. Each reaction requires 25 µl of the Taq OptiMix CLEAR. Simply add primers, template and water to a total reaction volume of 50 µl.

Taq OptiMix CLEAR offers several advantages: Reduced set up time and risk of contamination. Fewer reagent handling steps resulting in higher reproducibility. Optimized for increased specificity of DNA targets up to 4 kb.

Composition of the Taq OptiMix CLEAR (2x)

- Tris-HCl pH 8.5, (NH₄)₂SO₄, 3 mM MgCl₂, 0.2% Tween® 20
- 0.4 mM of each dNTP
- Ampliqon Taq DNA polymerase

Recommended Storage and Stability

Long term storage at -20 °C. Product expiry at -20 °C is stated on the label.

Quality Control

Taq DNA Polymerase is tested for contaminating activities, with no traces of endonuclease activity, nicking activity or exonuclease activity.

Protocol

This protocol serves as a guideline to ensure optimal PCR results when using Taq OptiMix CLEAR 2x Master Mix. Optimal reaction conditions such as incubation times, temperatures, and amount of template DNA may vary and must be determined individually.

1. Thaw Taq OptiMix CLEAR and primers. **It is important to thaw the solutions completely and mix thoroughly before use to avoid localized concentrations of salts.** Keep all components on ice.
2. Set up each reaction. Table 1 shows the reaction set up for a final volume of 50 µl. If desired, the reaction size may be scaled down. Use 12.5 µl Taq OptiMix CLEAR in a final volume of 25 µl.

Table 1. Reaction components (reaction mix and template DNA)

Component	Vol./reaction*	Final concentration*
Taq 2x Master Mix	25 µl	1x
25 mM MgCl ₂	0 µl (0-6 µl)	1.5 mM (1.5 – 4.5 mM)
Primer A (10 µM)	1 µl (0.5 – 5 µl)	0.2 µM (0.1 – 1.0 µM)
Primer B (10 µM)	1 µl (0.5 – 5 µl)	0.2 µM (0.1 – 1.0 µM)
PCR-grade H ₂ O	X µl	-
Template DNA	X µl	genomic DNA: 50 ng (10 – 500 ng) plasmid DNA: 0.5 ng (0.1 – 1 ng) bacterial DNA: 5 ng (1 – 10 ng)
TOTAL volume	50 µl	-

* Suggested starting conditions; theoretically used conditions in brackets

3. Mix the reaction mix thoroughly and dispense appropriate volumes into reaction tubes. Mix gently, e.g. by pipetting the reaction mix up and down a few times.
4. Add template DNA to the individual tubes containing the reaction mix.
5. Program the thermal cycler according to the manufacturer's instructions and recommendations in table 2. (DNA targets < 1kb) or table 3. (DNA targets 1 – 4 kb)
For maximum yield and specificity, temperatures and cycling times should be optimized for each new template target or primer pair.
6. Place the tubes in the thermal cycler and start the reaction.

Table 2. Three-step PCR program for targets < 1kb

Cycles	Duration of cycle	Temperature
1	2 – 5 minutes	95 °C
25 - 35	20 – 30 seconds ^a 20 – 40 seconds ^b 30 seconds ^c	95 °C 50 – 65 °C 72 °C
1	5 minutes ^d	72 °C

Table 3. Three-step PCR program for targets 1kb – 4kb

Cycles	Duration of cycle	Temperature
1	2 – 5 minutes	95 °C
25 - 35	20 – 30 seconds ^a 20 – 40 seconds ^b 60 – 300 seconds ^c	95 °C 50 – 65 °C 72 °C
1	5 minutes ^d	72 °C

^a. Denaturation step: This step is the first regular cycling event and consists of heating the reaction to 95 °C for 20 – 30 seconds. It causes melting of the DNA template by disrupting the hydrogen bonds

between complementary bases, yielding single-stranded DNA molecules.

- ^b Annealing step: The reaction temperature is lowered to 50 – 65 °C for 20 – 40 seconds allowing annealing of the primers to the single-stranded DNA template. Typically, the annealing temperature is about 3 – 5 °C below the T_m (melting temperature) of the primers used.
- ^c Extension/elongation step: Taq polymerase has its optimum activity temperature at 72 °C. At this step the DNA polymerase synthesizes a new DNA strand complementary to the DNA template strand. The extension time depends on the length of the DNA fragment to be amplified. As a rule of thumb, at its optimum temperature the DNA polymerase will polymerize a thousand bases per minute.
- ^d Final elongation: This single step is occasionally performed at a temperature of 72 °C for 5 minutes after the last PCR cycle to ensure that any remaining single-stranded DNA is fully extended.

Two-step PCR program

Fast 2-step PCR protocols are available using this link:
<https://ampliqon.com/en/pcr-technology/application-notes/>

Notes:

- The final $MgCl_2$ concentration of Taq OptiMix CLEAR is 1.5 mM. In some applications, more than 1.5 mM $MgCl_2$ is required for best results. Use 25 mM $MgCl_2$ to adjust the Mg^{2+} concentration according to table 4.

Table 4. Additional volume (µl) of $MgCl_2$ per 50 µl reaction:

Final $MgCl_2$ conc. in reaction (mM)	1.5	2.0	2.5	3.0	3.5	4.0	4.5
Volume of 25 mM $MgCl_2$	0	1	2	3	4	5	6

Related Products

Taq Master Mixes (500 x 50 µl reactions) *	Cat. No.
2x Master Mix, 1.5 mM $MgCl_2$ final concentration	A140303
2x Taq OptiMix CLEAR, 1.5 mM $MgCl_2$ final concentration	A370503
2x Master Mix RED, 1.5 mM $MgCl_2$ final concentration	A180303

TEMPase Hot Start Master Mixes (500 x 50 µl reactions) *	Cat. No.
2x Master Mix A**, 1.5 mM $MgCl_2$ final concentration	A230303
2x Master Mix A**BLUE, 1.5 mM $MgCl_2$ final concentration	A290403

*Master mixes available also in 1.1x variants as well as 2 mM $MgCl_2$ variants, **Mix A is Ammonium Buffer based, also available as Mix C based on Combination Buffer.

Special TEMPase Master Mixes (500 x 50 µl reactions)	Cat. No.
Multiplex 2x Master Mix, 3 mM $MgCl_2$ final concentration	A260303
GC TEMPase 2x Master Mix I – for GC-rich templates	A331703
GC TEMPase 2x Master Mix II – for GC-rich templates	A332703

Taq DNA Polymerase (500 units) *	Cat. No.
Taq DNA Polymerase 5 U/µl	A110003
• with 10x Ammonium Buffer	A111103

*Available in kits including one or two buffers (Ammonium Buffer, Standard Buffer or Combination Buffer). All kits include extra 25 mM $MgCl_2$

Hot Start DNA Polymerase (500 units) *	Cat. No.
TEMPase Hot Start DNA Polymerase, 5 U/µl	A220003
• with 10x Ammonium Buffer	A221103

*Available in kits including one or two buffers (Ammonium Buffer, Standard Buffer or Combination Buffer). All kits include extra 25 mM $MgCl_2$

Buffers for DNA polymerases *	Cat. No.
10x Ammonium Buffer, 3 x 1.5 ml	A301103
10x Standard Buffer, 3 x 1.5 ml	A302103
10x Combination Buffer, 3 x 1.5 ml	A303103
5x PCR Buffer RED, 6 x 1.5 ml **	A301810
PCR Grade Water, 6 x 5 ml	A360056

*Ammonium Buffer, Standard Buffer and Combination Buffer are also available as Mg^{2+} free buffers, detergent free buffers and Mg^{2+} and detergent free buffers.

**For direct gel loading and visualisation.

For Research Use Only. Not for use in diagnostics procedures.

Other product sizes, combinations and customized solutions are available. Please look at www.ampliqon.com or ask for our complete product list for PCR Enzymes. For customized solutions please contact us.

Made in Denmark

Issued 08/2021